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Instruktioner til SmartClip™ deaktiveringstang

3M Unitek
Orthodontic Products
2724 South Peck Road
Monrovia, CA 91016 USA

European Representative:
3M Unitek
Ohmstrasse 3
86899 Landsberg, Germany
+49/(0)8191/ 9474-0


Dentronix is a registered trademark of Dentronix, Inc.
Steril-Quik is a trademark. SmartClip is a trademark of 3M.

© 2006, 2009 3M. All rights reserved. REF 011-618-2 0904

SmartClip™ Appliance Wire Disengagement Instrument Instructions

Instructions For Use


English

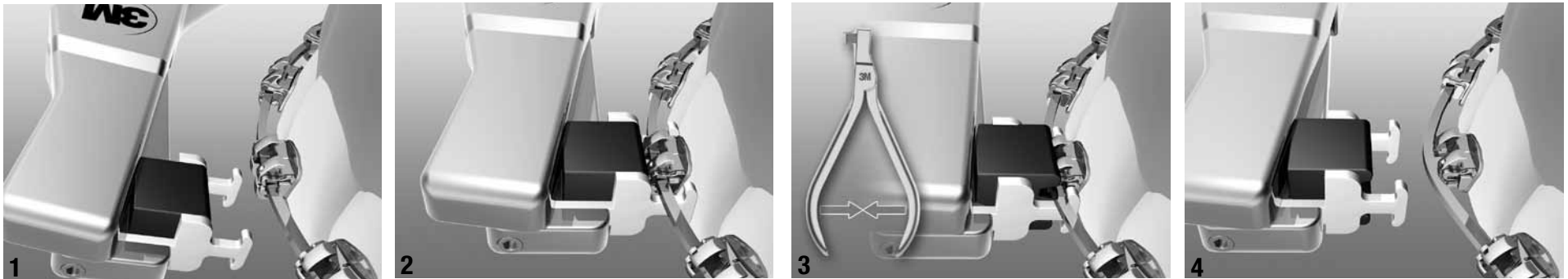Disengagement

1. Approach archwire with disengagement instrument from occlusal or gingival side. (Figure 1)
2. Place instrument hooks under archwire keeping the bracket between the instrument hooks. (Figure 2)
Make sure the instrument hooks are not under the bracket extensions holding the clips to avoid debonding or breakage.
The instrument hooks are to be mesial and distal of the bracket only.
3. Squeeze handles. (Figure 3)


Note: Repeat procedure for each bracket.

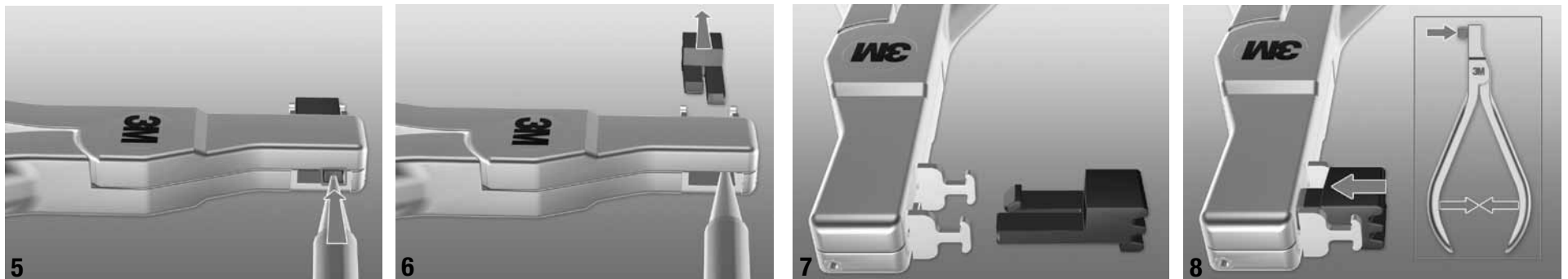Slide Replacement

1. Using a pointed object, approach tab.
2. Push tab. (Figure 5)
3. Remove slide. (Figure 6)
4. Obtain new slide. **Note:** Ensure proper orientation of slide. (Figure 7)
5. Push slide while holding handles in closed position. (Figure 8)

SmartClip™ Appliance Wire Disengagement Instrument Instructions

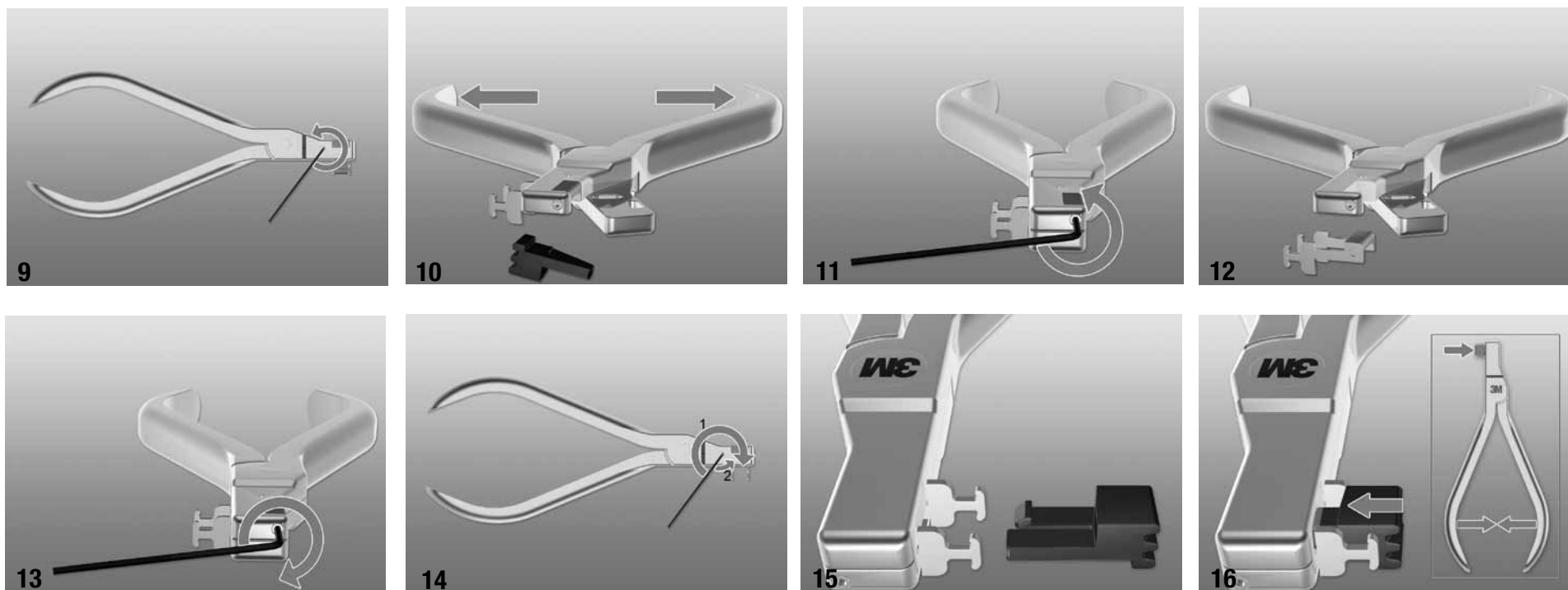Hook Replacement

1. Insert allen wrench* on first setscrew (Figure 9) and rotate counter-clockwise until handles open easily. (Figures 9 & 10)
2. Insert allen wrench on top setscrew and rotate counter-clockwise until hooks are free. (Figures 11 & 12)
3. Obtain new hooks. **Note:** Ensure proper orientation of hooks. (Figure 12)
4. Tighten top setscrew until first resistance is felt. (Figure 13)
5. Close plier handles. Insert allen wrench and tighten first setscrew clockwise until first resistance is felt and then rotate counter-clockwise 1/4 turn. (Figure 14)
6. Obtain new slide. **Note:** Ensure proper orientation of slide. (Figure 15)
7. Push slide while holding handles in closed position. (Figure 16)

*Included with Hook Insert Replacement Kit (P/N 804-162)

SmartClip™ Appliance Wire Disengagement Instrument Instructions

WARNINGS:

Use instruments for indicated use only. Failure to do so may result in the breakage of the instrument and injury to patient or orthodontic professional. Safety glasses recommended for patient and orthodontic staff when using these instruments.

Stainless steel contains nickel and chromium. If an allergic reaction occurs, direct patient to consult a physician.

Disinfection Information:

Please refer to literature 011-650, “Reprocessing Instructions for Use”.

Warranty Information

3M Unitek warrants that this instrument will be free from defects in materials, or workmanship for a period of 1 year from the date of purchase. Failure to follow recommended use, maintenance and sterilization procedures will void all warranty claims.

Recommended maintenance includes cutter tip/blade re-sharpening (if applicable) and reconditioning (cleaning, polishing, adjusting screw as necessary) every six months. 3M Unitek’s sole obligation and customer’s sole remedy in the event of a claimed defect shall be limited to, at 3M Unitek’s option, repair, replacement of merchandise, or refund of the purchase price. ALL IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE OR USE ARE DISCLAIMED.

Limitation of Liability

Except where prohibited by law, 3M Unitek will not be liable for any loss or damages arising from this Product, whether direct, indirect, special, incidental or consequential, regardless of the theory asserted, including warranty, contract, negligence or strict liability. This limitation does not apply to third party personal injury claims.

All instruments should be sent to the 3M Unitek Authorized Service Center at the address below for evaluation to determine credit, replacement, or repair. The Authorized Service Center will determine, upon receipt of instrument, if the instrument is to be covered by warranty. Materials, workmanship, and corrosion resistance are the only items covered by this warranty. Misuse or abuse will void warranty. Customer will be contacted regarding repair charges for out of warranty instruments.

3M Unitek – Authorized Service Center

Precision Plier Service

14175 Telephone Avenue, Suite D

Chino Hills, CA 91710-5762

Phone Number: 1-800-229-9060 USA, Canada. 1-909-590-2085.

References

ISO 17664:2004(E), Sterilization of medical devices—Information to be provided by the manufacturer for the processing of resterilizable medical devices.

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